

Provox Life™ Stability Adhesive Technical Data

Provox Life™ Stability Adhesive

Technical Data

Provox Life™ Standard Adhesive







Provox Life™ Sensitive Adhesive



Provox Life™ Night Adhesive



Provox Life™ Stability Adhesive

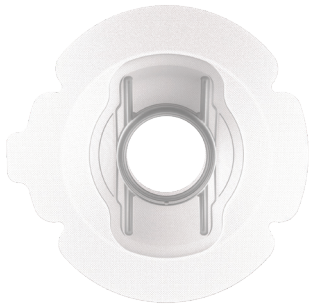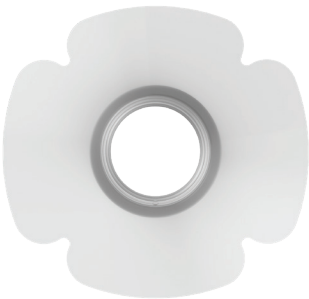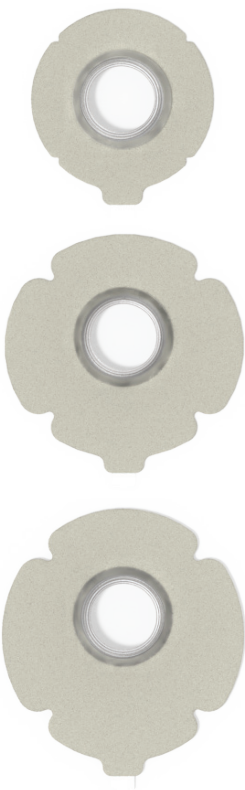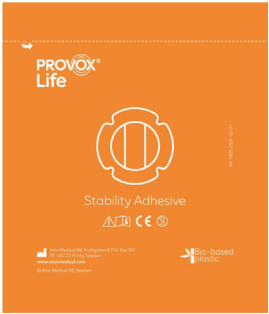

Provox Life™ Stability Adhesive

Technical Data

Carrier:

Polyethylene with siliconised surface. Perforated.

Adhesive:

Medical Grade Acrylic

Liner:

Polypropylene. Three-piece.

Connector:

SecureFit. Diameter 23mm.

Firm, concave base with stabilising bars for deeper stomas and neck support during (handsfree) TE speech.

Application:

The pressure sensitive medical grade acrylic material reaches maximum adhesion after about one hour and does not benefit from pre-warming prior to application.

Removal:

Reduced moisture transport rate (MTR) may lead to higher levels of moisture trapping underneath the adhesive which helps to reduce stickiness for adhesive removal. Use Provox Adhesive Remover wipe on top of perforated carrier to saturate material and loosen adhesive prior to removal. Then use adhesive remover wipe from the side to gently remove adhesive – starting from finger lift tab.

Technical data versus legacy device Provox Stabilibase adhesive:

The design of the firm, concave base and connector has remained the same. Materials same as Provox Life Standard.

Sizes:

Available in one size only.

At a glance:

Provox Life Standard Adhesive Provox Life Stability Adhesive	Material: Medical grade acrylic
Flexibility	High
Permeability	Moderate
Absorption	None
Adhesion	High tack
Risk of skin stripping	High, if not removed correctly